

AERO RUNNER® TREADMILL 2.0

The new generation of Anti-Gravity Treadmill—a comprehensive upgrade—leads gait rehabilitation.

ServoStep

Wider adjustable speed range:
-10 km/h to 30 km/h; (GA300/GA300S Only)
Minimum starting speed: 0.1 km/h;
Adjustment precision: 0.1 km/h.

Longer armrests and wider
armrest spacing



Lower step-in height for easier
patient access

Heart Rate Controlled Running (Optional)

Using a heart rate strap,
the treadmill automatically
changes speed to keep
your heart rate at a set
value, enabling precise
rehab training.




Electric lift armrests for easy
height adjustment
(GA100S/GA200S/GA300S Only)

Perturbation mode (GA100S/GA200S/GA300S Only)


Randomly generated speed
fluctuations simulate complex
real-world walking environments
to train patients' balance regula-
tion ability.

MORE DETAILED REHABILITATION DATA, MORE DIVERSE TRAINING MODES


Center of Pressure (COP) measure

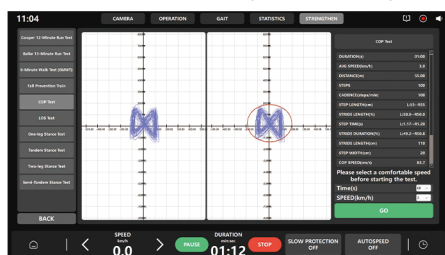

Dynamic and Static Balance Training

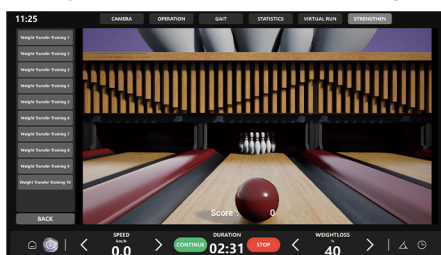

Surround Sound Training

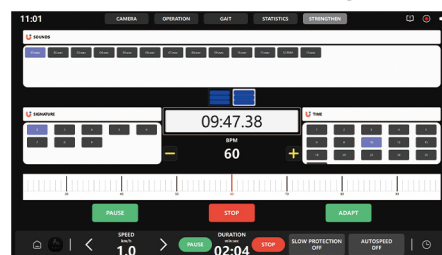

UNLOCK NEW POSSIBILITIES



Perturbation Training Board



Situational Interaction System



Supports VR Headset
Interactive Scenario Mode

WHY CHOOSE AERO RUNNER?

Innovative Efficiency for Exceptional Value

- **Complete In-House Development**

Ensures top-quality products tailored to your rehabilitation needs and continuous improvement and update of software and hardware.

- **Focus on the Long-term Interests of Users**

You can get the best return on investment and experience throughout the user's life cycle.

Unparalleled Service and Support

- **Simplified Setup**

Shipment in whole units for easier installation.

- **Industry-Leading Warranty**

2 years free.

Up to 10 years for ultimate peace of mind. (Optional)

- **Personalized Touch**

Customized language.

	GA300	GA200	GA100
	The GA100S, GA200S, GA300S models offer exclusive features and accessory upgrades		
	Report Printing ★ Perturbation mode ★ Armrests Motorized adjustment Touchscreen Size: 32 inches		
★ User Weight Range	36 kg – 185 kg		
★ Weight Support Adjustment	0–100% range with 1% adjustment steps		
★ Speed	-10 km/h to 30 km/h	-10 km/h to 25 km/h	-10 km/h to 20 km/h
Starting Speed	0.1 km/h		
Speed regulation unit	0.1 km/h		
★ Incline	Slope can be adjusted in 19 levels from 0 to 18% (0–10°)		
Live Video Monitoring	3-way simultaneous monitoring (Left, Right, Front)		
Infrared Gait Sensor	Standard		
Gait Display	Standard		
Airbag Material	Environmentally friendly and odorless		
★ Center of Gravity Transfer Balance Games	10 types		
★ Heart Rate Monitoring and Heart Rate Controlled Running	Optional		
Warranty	Up to 10 years (Optional)		
Situational Interaction System	4 different scenarios (Optional)		
★ Touchscreen Size	23.8 inches		
Armrest Adjustment	Manual adjustment		

